

Non-target screening of new compounds in Arctic apex predators

Mikael Harju, Laura Røhler, Pawel Rostkowski

NILU

Heli Routti, Jon Aars, Christian Lydersen

Kit M. Kovacs, Geir Wing Gabrielsen

Norwegian Polar institute




Aim of the study

- The effects of Contaminants on polar bear energetics – using *in vitro* exposure studies on the peroxisome proliferator-activated receptors (PPARs)

Binding \leftrightarrow Activation

- Extracts as well as synthetic mixtures was compared
- Screening of potential new pollutants

Project outline

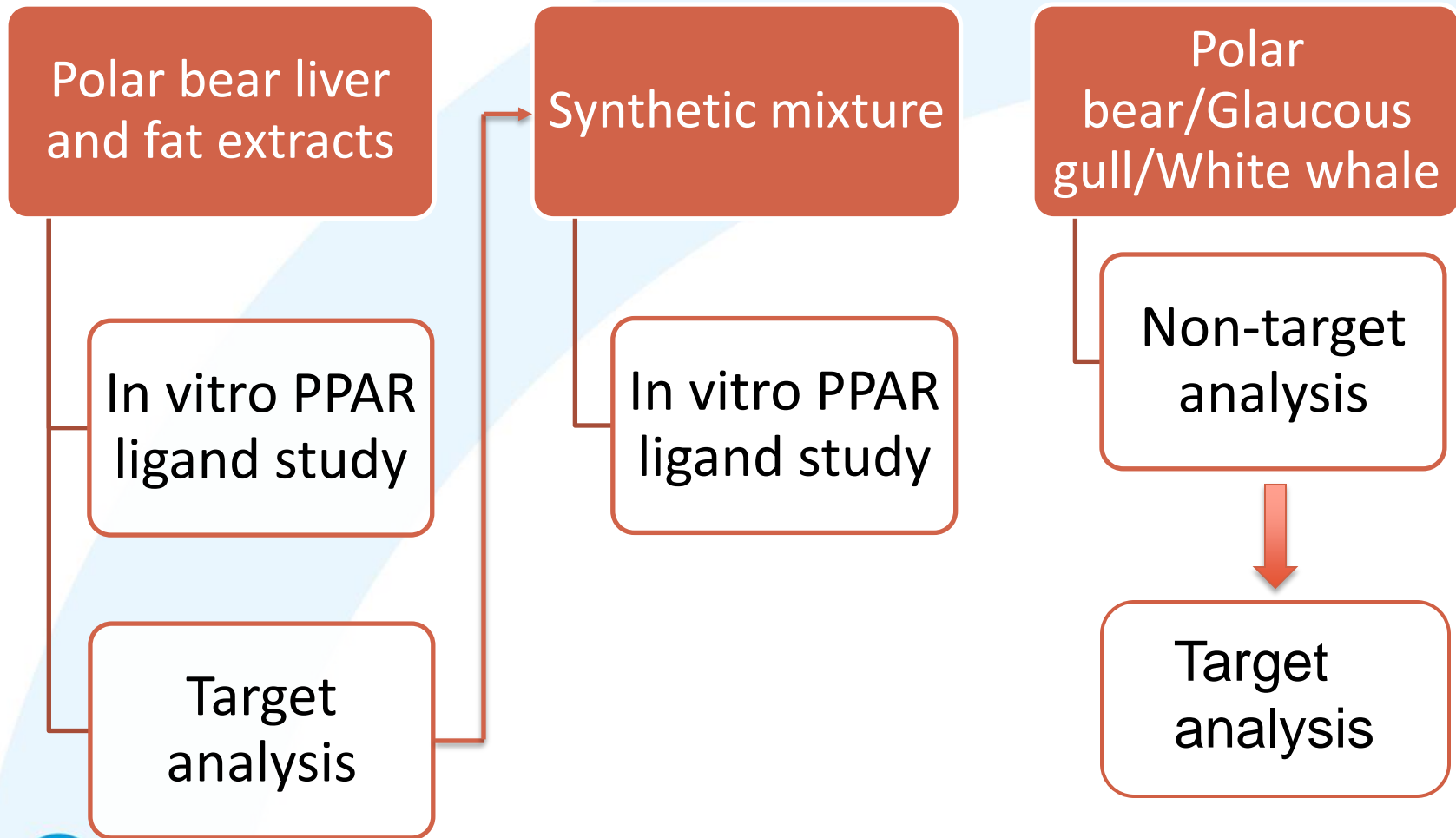

Extraction

68g Polarbear liver
12g blubber
Control/blank + QC

Homogenization
Dry Na_2SO_4

Extracting
Lipids and POPs


Open-column Extraction



Cleanup

Reverse Osmosis

Thin layflat polyethylene tubing
10 \AA pores

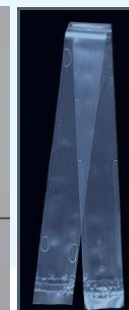

Removes:
88% lipids liver
>98% in blubber

GPC/SEC
Size exclusion
chromatography

Chromatography
On Florisil

F1

F2


F3

F1: Neutrals PCBs,
PBDEs, pest etc
F2: MeSO_2 -PCBs,
 MeSO_2 -DDE
F3: OH-PCBs, PCP,
OH-PBDEs

Conc factor 100 to 1

Target screening

PCB composition

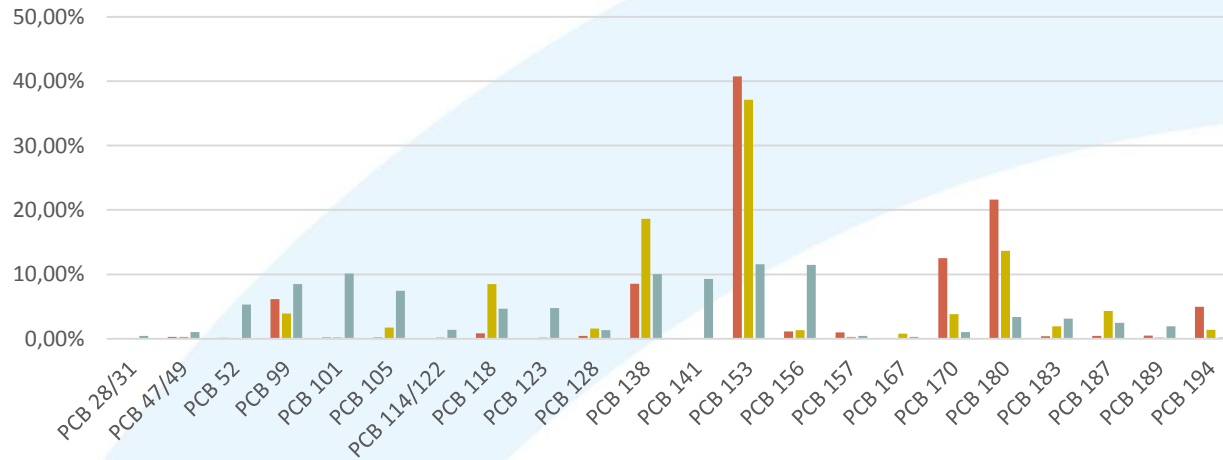

Biotransformation :

■ Polar bear > ■ Glaucous gull > ■ White whale




Non-target analyses

- Fat samples from 3 polar bears, 3 glaucous gulls, 3 White whales
- Svalbard 2013-14
- Extraction: Homogenized in Na_2SO_4 followed by extraction with acetone/cyclohexane
- Clean-up: semipermeable membrane, gel permeation chromatography and fractioning on Florisil to three according to polarity of the compounds
- Analyses: Leco Pegasus 4D, 2-D gas chromatografi coupled to time of flight mass spectrometry
- Data processing: Leco ChromaTOF software

Peak find → Peak deconvolution → matching to spectra database (NIST)
→ scripting

White whale

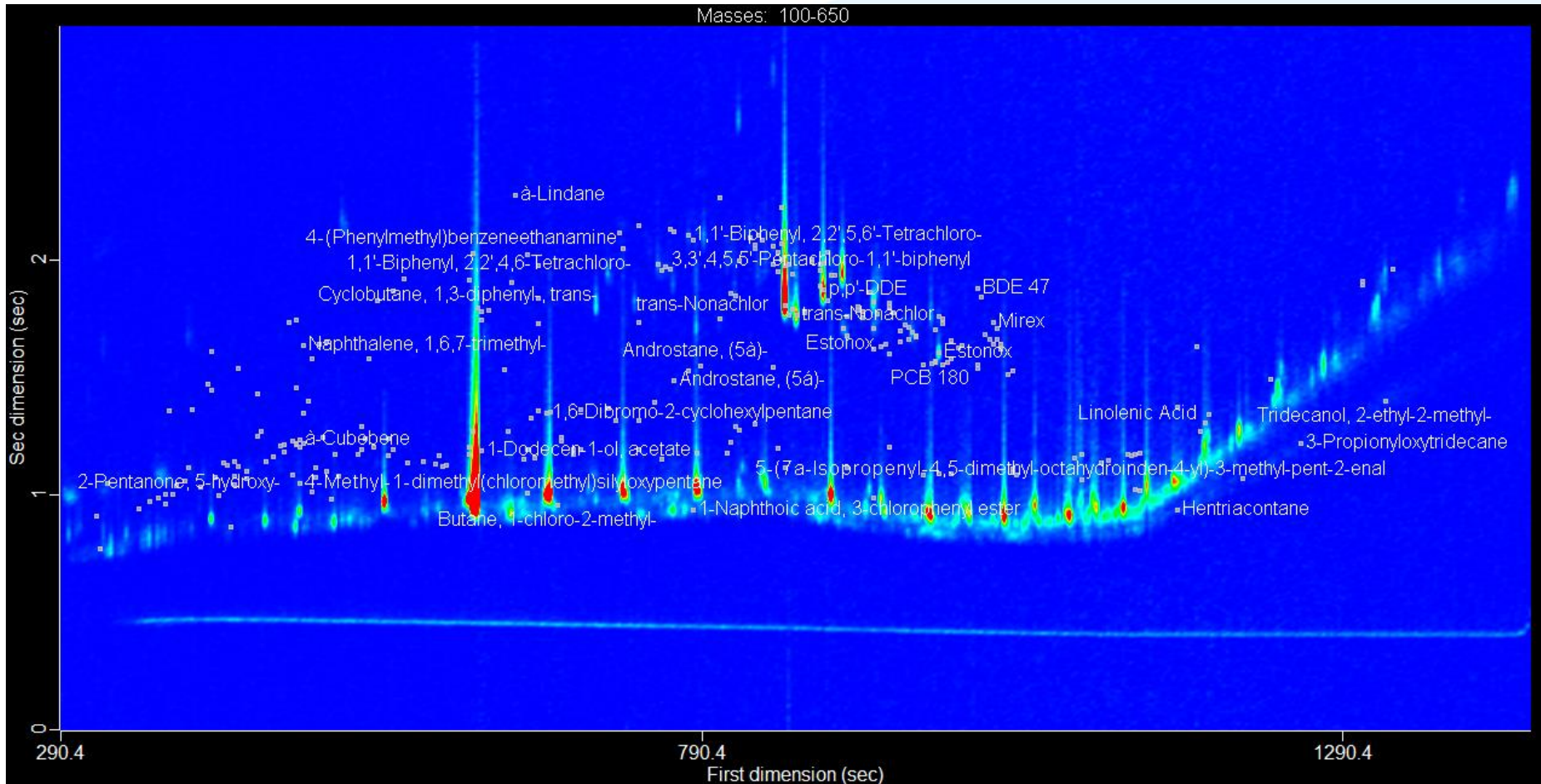

Results of new compounds by non-target analyses

Hit	CAS	Polar bear	Glaucous gull	White whale
Octachlorstyrene	29082-74-4	x	x	
Heptachlor-endo-epoxide (isomer A)	28044-83-9	x	x	
2,3-Dichlorobenzo[b]thiophene	5323-97-7	x	x	
1-Piperidinecarbodithioic acid, 2,3,5,6-tetrachloro-4-pyridinyl ester		x		
(2,3-Diphenylcyclopropyl)methyl phenyl sulfoxide, trans-	131758-71-9	x		x
9-Allylanthracene	23707-65-5			x
l-Proline, n-pentafluoropropionyl-, heptyl ester				x
2-Chloro-1,1-bis(4-chlorophenyl)-ethene (DDMU)	1022-22-6		x	x
7-Methylbenz(a)anthracene	2541-69-7			x
Tris(3-chlorophenyl)phosphine	29949-85-7			x
6-Methoxy-2,2',4,4'-tetra (BC-3)	102739-99-1			x
Naphthalene, 1,2,3,4-tetrahydro-1-phenyl-				x
Bipyrrol Q1	428442-17-5			x
1,3-Cyclopentadiene, 1,2,3,4-tetrachloro-	695-77-2		x	
1-Silacyclopentadiene, 3-(diethylboryl)-4-ethyl-1,1-dimethyl-2-(3-methylbutyl)-5-trimethylstannyl-	153331-61-4		x	
Endosulfan	115-29-7		x	
2,3,4,5,6-Pentachlorophenyl acetat	1441-02-7		x	
2,7:3,6-Dimethanonaphth[2,3-b]oxiren-8-ol, 3,4,5,6,9,9-hexachloro-1a,2,2a,3,6,6a,7,7a-octahydro-, acetate, stereoisomer	55649-85-9		x	
1,4:5,8-Dimethanonaphthalene-2,3-diol, 5,6,7,8,9,9-hexachloro-1,2,3,4,4a,5,8a-octahydro-, diacetate, (1à,2à,3á,4à,4aá,5á,8á,8aá)-	34408-22-5		x	
Chlordecone	55570-85-9		x	

- **Qualitative**
- **No information of concentrations**
- **Structures not confirmed by standards**

Target analysis

- Fat samples from 10 polar bears, 8 glaucous gulls, 9 white whales
- Extraction and cleanup similar as non-target analysis (no semipermeable membrane)– full set of internal standards
- Analysis on GC-MS




Target analyses of compounds found in non-target screening:

Concentrations (ng/g ww)	Polar bear	Glaucous gull	White whale
Octachlorstyrene	8.7	21.6	3.8
Heptachlor-endo-epoxide	< 0.9	< 0.9	< 0.9
2,3-Dichlorobenzo[b]thiophene	< 0.1	< 0.1	< 0.1
1-Piperidinecarbodithioic acid, 2,3,5,6-tetrachloro-4-pyridinyl ester			
(2,3-Diphenylcyclopropyl)methyl phenyl sulfoxide, trans-			
9-Allylanthracene			
l-Proline, n-pentafluoropropionyl-, heptyl ester			
2-Chloro-1,1-bis(4-chlorophenyl)-ethene (DDMU)	2.5	54.3	3.5
7-Methylbenz(a)anthracene	< 0.1	< 0.1	< 0.1
Tris(3-chlorophenyl)phosphine	< 0.3	< 0.3	< 0.3
6-Methoxy-2,2',4,4'-tetra (BC-3)			
Naphthalene, 1,2,3,4-tetrahydro-1-phenyl-	< 2.5	< 2.5	< 2.5
Bipyrrol Q1			
1,3-Cyclopentadiene, 1,2,3,4-tetrachloro-	< 0.3	< 0.3	< 0.3
1-Silacyclopentadiene, 3-(diethylboryl)-4-ethyl-1,1-dimethyl-2-(3-methylbutyl)-5-trimethylstannyl-			
Endosulfan	Will be analyzed		
2,3,4,5,6-Pentachlorophenyl acetat			
2,7:3,6-Dimethanonaphth[2,3-b]oxiren-8-ol, 3,4,5,6,9,9-hexachloro-1a,2,2a,3,6,6a,7,7a-octahydro-, acetate, stereoisomer			
1,4:5,8-Dimethanonaphthalene-2,3-diol, 5,6,7,8,9,9-hexachloro-1,2,3,4,4a,5,8,8a-octahydro-, diacetate, (1à,2à,3á,4à,4aá,5á,8á,8aá)-			
Chlordecone	< 0.1-2.5	24.7	1.7

- **New extractions**
- **Quantification using reference standards**
- **Retention time confirmed using authentic standards**
- **Quantitative results**

Target analyses of compounds found in non-target screening:

- Octachlorostyrene:
 - Industrial by-product, not regulated
 - Not new, reported in polar bears already in 2000
 - Highly toxic
- DDMU:
 - Metabolite of DDE
 - Toxic
- Chlordecone/Kepone
 - Pesticide, included in Stockholm Convention in 2009
 - Not measured in the Arctic (Muir & deWit, 2010)
 - Toxic

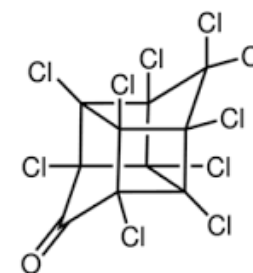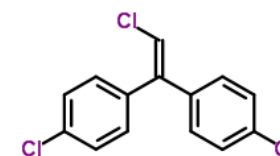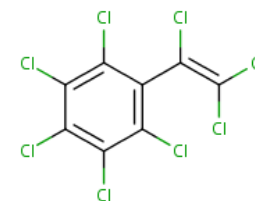

The way forward

- Suspect screening, good alternative or complement to non-target screening
- Accurate high resolution MS for additional confirmation - standards are difficult and expensive to get
- Confirmation of compounds – Retention time, accurate mass, mass fragmentation spectra, isotope ratios
- Negative and positive chemical ionization
→ Selectivity and molecular fragments

For more information, www.nilu.no

